

PIXIS 100/400 CCD CAMERA

KEY FEATURES

- 1340 x 100 format sensor for spectroscopy
- Large 20 μm pixels to maximize sensitivity
- eXcelon® technology for sensitivity across UV and NIR
- Thermoelectric cooling down to 90 °C for low dark current
- Acquisition rates of up to 2 MHz
- Permanent vacuum seal
- Dual amplifiers with software selectable gains
- Optional UV phosphor coatings
- Optional AR coating and wedge windows are available

TYPICAL APPLICATIONS

- Raman spectroscopy
- Single-molecule spectroscopy
- Absorbance spectroscopy
- Emission spectroscopy
- Fluorescence spectroscopy
- Reflectance spectroscopy

RELIABILITY

- One-year warranty
- Lifetime vacuum guarantee

High-Performance, Low-Noise CCDs

The PIXIS camera family consists of fully integrated, low noise CCD cameras with deep cooling for quantitative scientific imaging and spectroscopy applications from UV to NIR. High quantum efficiency and ultra-low-noise electronics make PIXIS cameras highly sensitive and ideal for demanding low-light applications.

These models include the PIXIS 100 and 400 CCDs, featuring a back-illuminated 1340 x 100 or 1340 x 400 sensor format, respectively, and large 20 μm pixels optimized for sensitive spectroscopy research. The PIXIS 100 and 400 offer thermoelectric cooling down to -80 °C, an all-metal hermetically sealed design (with a lifetime vacuum guarantee), and exclusive eXcelon® technology for increased sensitivity from UV to NIR with reduced etaloning.




PIXIS 100/400 SPECIFICATIONS

SPECIFICATIONS	PIXIS 100	PIXIS 400
CCD Sensor Format	B: Back-illuminated (highest QE in visible, subject to NIR etaloning, optional UV coating) BR: Back-illuminated, deep-depleted (high QE in NIR with minimal etaloning, optional UV coating) BX: Back-illuminated, with eXcelon (high QE from UV-NIR with low etaloning) BRX: Back-illuminated, deep-depleted with eXcelon (highest QE from UV-NIR with negligible etaloning)	
Active Array Size	1340 x 100 pixels	1340 x 400 pixels
Pixel Area	20 μm x 20 μm (400 μm^2)	
Sensor Area	26.8 mm (W) x 2.0 mm (H) (optically centered)	26.8 mm (W) x 8.0 mm (H) (optically centered)
Spectral Response	250 – 1100 nm	
Spectrometric Well Capacity	300,000 e ⁻ (High Sensitivity); 1,000,000 e ⁻ (High Capacity)	
ADC Rates	2 MHz, 100 kHz	
Analog Gain (typical)	1, 2, 4 e ⁻ /ADU (High Sensitivity); 4, 8, 16 e ⁻ /ADU (High Capacity) (available at all ADC rates)	
Non-linearity	< 1.0% @ 100 kHz	
I/O Signals	Two MCX to BNC connectors for programmable frame readout, shutter, trigger in	
Typical System Readout Modes	3.0 e ⁻ @ 100 kHz; 11.0 e ⁻ @ 2 MHz	
Typical Dark Current	BRX/BR: 0.03 e ⁻ /p/s @ -80 °C BX/B: 0.001 e ⁻ /p/s @ -80 °C	
Cooling Options	Thermoelectric air or liquid cooling (optional CoolCUBE II coolant circulator available) -70 °C	
Thermostating Precision	± 0.05 °C across temperature range	
Vertical Shift Times	< 15 μsec /row (programmable)	
Output Interface	USB 2.0	
Digitization	16-bit	
Operating Temperature Range	5 °C to 30 °C (non-condensing)	
Camera Weight	5 lbs, 2.27 kg	
Camera Dimensions	163 x 118 x 114 mm (L x W x H)	
Certification	CE	


PIXIS 100/400 SPEED TABLES

SPECIFICATIONS	PIXIS 100	PIXIS 400
@100 kHz Full Vertical Bin	70 fps	60 fps
@2 MHz Full Vertical Bin	750 fps	315 fps
@2 MHz 0.4 mm high	1300 fps	1300 fps


PIXIS 100/400 QE CURVES

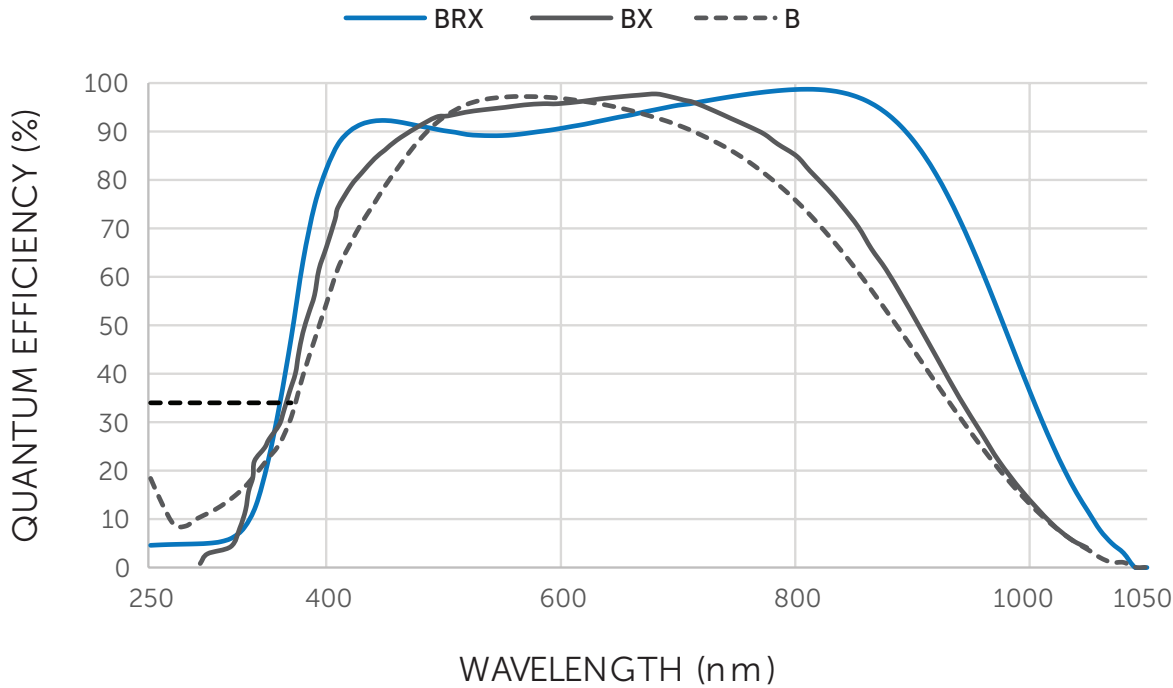

NIR FRINGE COMPARISON

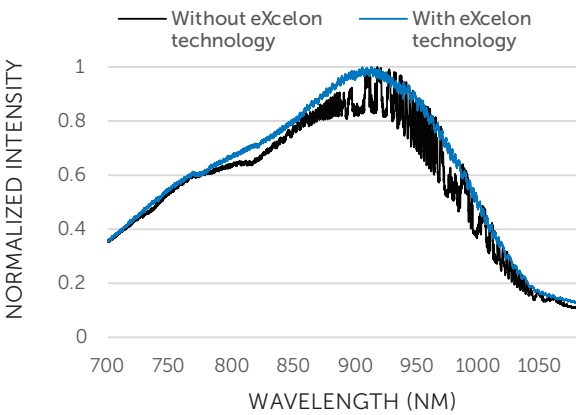


EXCELON QE IMPROVEMENT

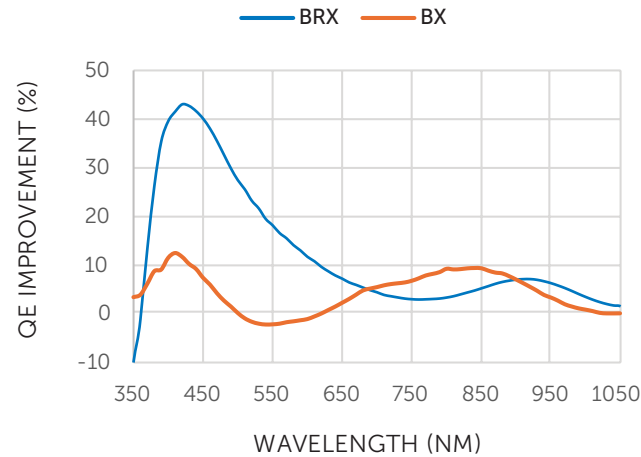


PIXIS 100/400 ACCESSORIES

Accessories (Air Cooled)	Accessories (Liquid Cooled)	Accessories (Optional)
MCX to BNC cables (x2)	MCX to BNC cables (x2)	WinX application software
USB 2.0 interface cable	USB 2.0 interface cable	LightField® application software
Certificate of Performance	Certificate of Performance	SITK™ for LabVIEW
Power supply (120V/240V)/cables	Power supply (120V/240V)/cables	Internal 25 mm or 45 mm shutter
	Liquid circulator with hoses	Fiber optic extender kit

PIXIS 100/400 DIMENSIONAL OUTLINES — SPECTROMETER MOUNT (UNITS: INCHES [MM])

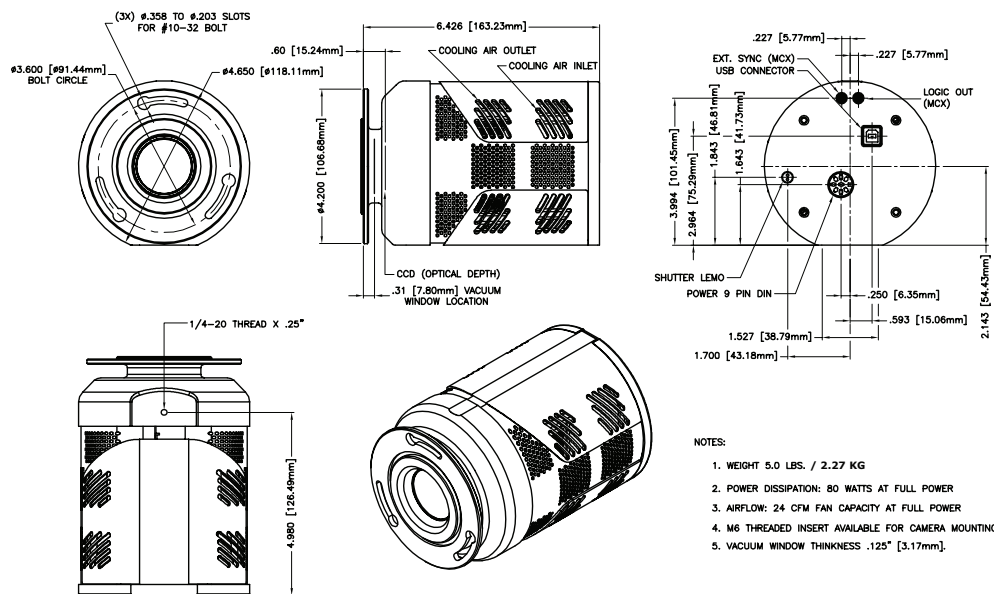

PIXIS 100/400 DIMENSIONAL OUTLINES — SPECTROMETER MOUNT WITH SHUTTER (UNITS: INCHES [MM])

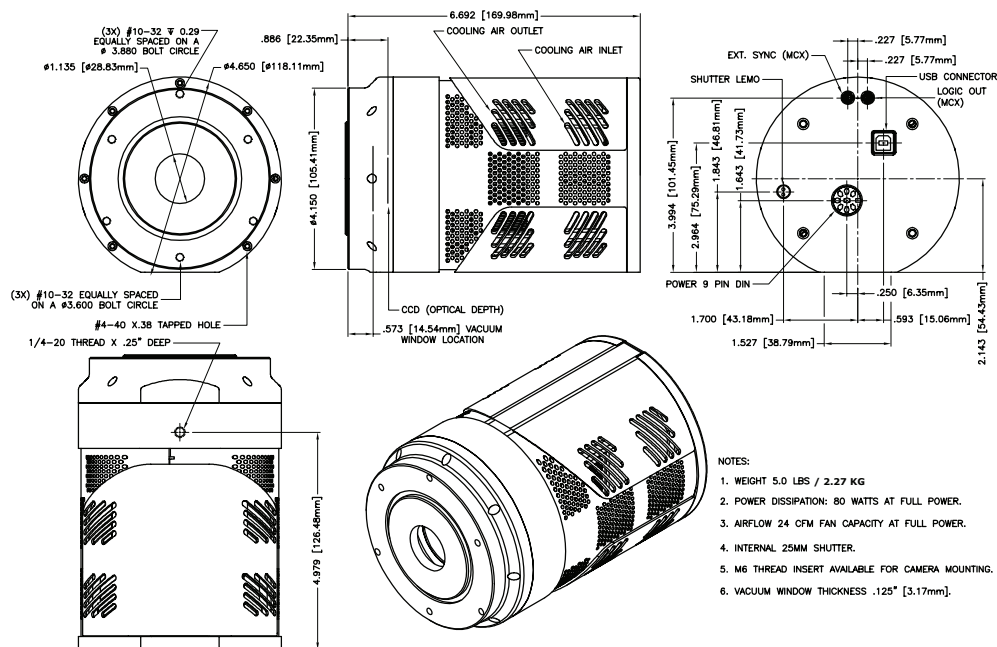

FOR MORE INFORMATION CONTACT:

AMERICAS	Trenton, New Jersey, USA +1 609-587-9797 pi.info@teledyne.com
EUROPE	France +33 (0)1 70 38 19 00 evr@teledyne.com
	Germany +49 (0)89 660 779 30 pi.germany@teledyne.com
	United Kingdom +44 (0) 7810 835 719 pi.info@teledyne.com
ASIA PACIFIC	China +86 157 2153 5343 pi.info.china@teledyne.com
	Japan +81 3 6709 0631 pi.nippon@teledyne.com

Specifications in this datasheet are subject to change.
Refer to the Teledyne Princeton Instruments website for
most current specifications.

Copyright © 2025 Teledyne Princeton Instruments, Inc.
All rights reserved.

Revision Date: 2025 01 24